

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 1 of 24

Contents

1	Introduction	3
1.1	Study Specifications and Acceptance Criteria	3
2	Materials and Equipment Description	3
3	Exploratory Data Analysis	3
4	Results	4
4.1	Within-Laboratory Precision Profile	4
4.2	LoQ Results for Simple Linear Regression with Variables Transformations . . .	5
4.3	LoQ Results for non-linear models	6
4.4	Results Summary	8
5	Software	9
6	Operator and Data Input Reviewer	10
7	Summary Raw Data	11

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 2 of 24

List of Figures

1	Precision profile plot with fitted linear regression model with lowest RMSE for lot N001.	6
2	Precision profile plot with fitted non-linear model with lowest AIC for lot N001.	8

List of Tables

1	LoQ Specification and Acceptance Criteria	3
2	Low Level Samples Exploratory Data Analysis Summary for lot N001.	4
3	Within-laboratory precision for low level samples for lot N001.	5
4	List of the models for multiple linear transformations for lot N001	6
5	List of the non-linear models fitted to the precision profile for lot N001.	7
6	Summary of the Raw Data for lot N001	11

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 3 of 24

1 Introduction

This is the report of the determination of the Limit of Quantitation (LoQ) for FVIII_Low_Test assay with lot N001 in the ModTOP platform. The report has been automatically generated. Data analysis has been executed by Joey. The LoQ study was performed following the approved guideline CLSI EP17-A2 Evaluation of Detection Capability for Clinical Laboratory Measurement Procedures, 2nd Edition.

A Functional Sensitivity approach determined the LoQ. This approach characterizes the precision profile in the lower end of the measuring interval and defines the LoQ as the lowest measurand concentration that achieves a desired within-laboratory precision or precision goal. The experiment tested 11 low level samples for 5 days, 1 runs per day and 10 replicates per run.

1.1 Study Specifications and Acceptance Criteria

LoQ has acceptance criteria and specification. LoQ results were acceptable if they met the acceptance criteria shown in Table 1.

Table 1. LoQ Specification and Acceptance Criteria

Specification (activity)	Acceptance Criteria (precision goal as %CV)
≤1.00	20.00

Abbreviations: CV, Coefficient of Variation.

2 Materials and Equipment Description

The materials and equipment used in the present LoQ study are described in the corresponding section of DRC-248-01.

3 Exploratory Data Analysis

The exploratory data analysis helps to have a first impression about the data.

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 4 of 24

Table 2. Low Level Samples Exploratory Data Analysis Summary for lot N001.

Sample	N	Min (activity)	Max (activity)
DP	50	0.00	0.20
A	50	0.30	0.80
B	50	0.50	1.10
C	50	0.70	1.40
D	50	1.10	1.90
E	50	1.30	2.20
F	50	1.90	2.80
G	50	2.20	3.00
H	50	2.60	3.40
I	50	3.00	4.50
J	50	3.60	4.70

Abbreviations: N, Sample Size; Min, Minimum; Max, Maximum.

4 Results

Functional Sensitivity approach characterized the precision profile of the low level samples data to determine the LoQ as the lowest measurand concentration that achieves a within-laboratory precision goal.

4.1 Within-Laboratory Precision Profile

A Variance Components Analysis (VCA) determined the within-laboratory precision for each sample via restricted maximum likelihood (REML). This method is robust to balanced and unbalanced study designs since the final data set might be unbalanced due to missing values or excluded results as statistical outliers.

Table 3 contains the VCA results summary with the mean, the sample size (N), the within-laboratory standard deviation and the within-laboratory coefficient of variation for low level samples data.

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 5 of 24

Table 3. Within-laboratory precision for low level samples for lot N001.

Sample	N	Mean (activity)	SD _{WL} (activity)	CV _{WL} (%)
DP	50	0.13	0.06	47.59
A	50	0.52	0.13	25.82
B	50	0.82	0.12	14.04
C	50	1.12	0.15	13.07
D	50	1.42	0.18	12.41
E	50	1.88	0.21	10.92
F	50	2.32	0.20	8.59
G	50	2.59	0.17	6.60
H	50	3.05	0.20	6.60
I	50	3.82	0.30	7.85
J	50	4.13	0.24	5.88

Abbreviations: N, Sample Size; SD_{WL}, Within-Laboratory Standard Deviation; CV_{WL}, Within-Laboratory Coefficient of Variation.

4.2 LoQ Results for Simple Linear Regression with Variables Transformations

This analysis fitted simple linear regression to the precision profile of within-laboratory coefficient of variation (CV_{WL}) versus mean concentration of low level samples (Figure 1), assessing multiple transformations of both variables.

Table 4 summarizes the assessed linear transformations with their equations, model parameters, goodness of fit parameters of root mean squared error (RMSE) and coefficient of determination (R²) and the estimated LoQ as the mean concentration at the CV_{WL} goal, sorted from lowest to greatest RMSE. Estimated LoQ at CV_{WL} goal (last column, Table 4) as N/A indicates a negative estimated LoQ or that the LoQ cannot be estimated. Figure 1 shows the fitted line of the lowest RMSE model along with its estimated LoQ.

R² utilizes the transformed Y scale (CV_{WL}), while RMSE the original Y scale. Selecting among models with different transformations of Y makes the RMSE the most meaningful goodness of fit parameter. The lower RMSE, the better fitting.

The linear precision profile model with the lowest RMSE interpolated a LoQ of 0.58 activity.

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 6 of 24

Table 4. List of the models for multiple linear transformations for lot N001

Equation	Transformation	β_0	β_1	R^2	RMSE	LoQ at CV_{WL} goal (activity)
$CV_{WL} = \beta_0 Mean^{\beta_1}$	$log(CV_{WL}), log(Mean)$	2.665	-0.6029	0.958	1.847	0.58
$CV_{WL} = \beta_0 + \beta_1 / Mean$	$1 / Mean$	7.304	5.32	0.936	2.986	0.42
$CV_{WL} = \beta_0 + \beta_1 \cdot log(Mean)$	$log(Mean)$	18.49	-11.43	0.913	3.462	0.88
$CV_{WL} = \beta_0 + \beta_1 Mean^{1/2}$	$Mean^{1/2}$	40.8	-20	0.727	6.140	1.08
$CV_{WL} = exp(\beta_0 + \beta_1 Mean)$	$log(CV_{WL})$	3.284	-0.4186	0.765	7.230	0.69
$CV_{WL} = (\beta_0 + \beta_1 Mean)^2$	$CV_{WL}^{1/2}$	5.204	-0.8148	0.658	7.386	0.90
$CV_{WL} = \beta_0 + \beta_1 Mean$	None	28.08	-6.854	0.544	7.940	1.18
$CV_{WL} = \beta_0 + \beta_1 Mean^2$	$Mean^2$	21.15	-1.203	0.330	9.624	0.98

Abbreviations: N/A, Not Applicable; R^2 , Coefficient of Determination; RMSE, Root Mean Squared Error; CV, Coefficient of Variation; WL, Within-Laboratory precision.

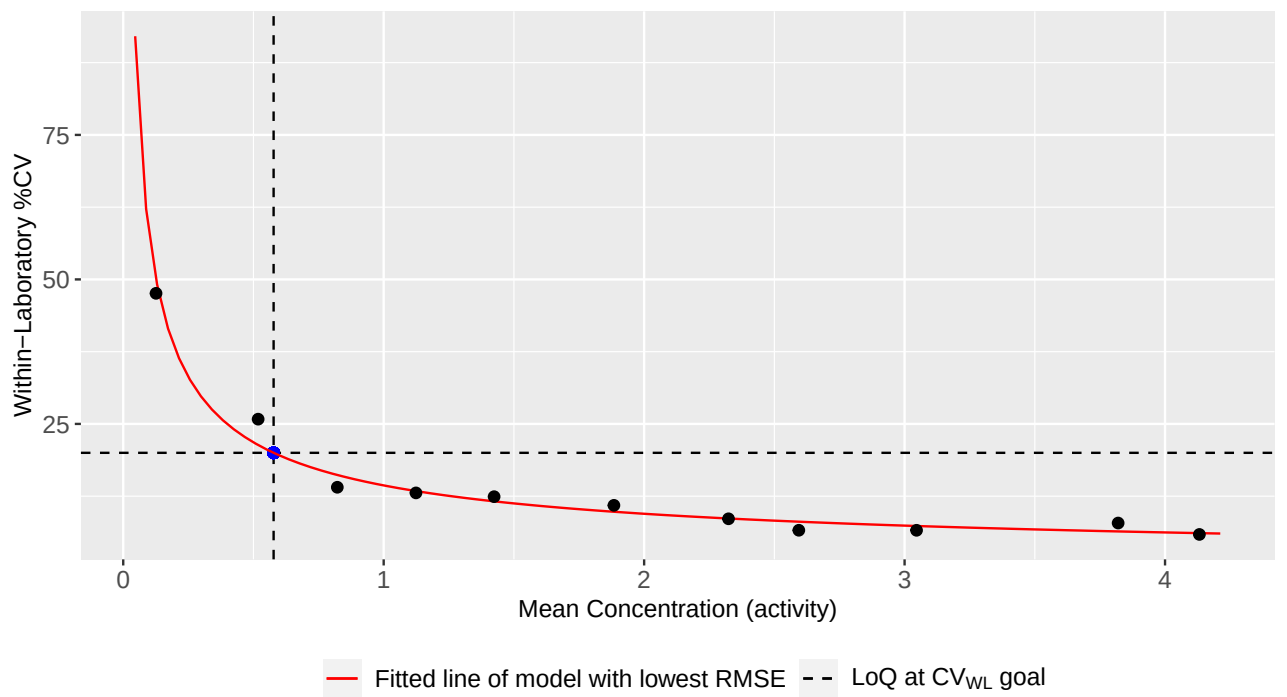


Figure 1. Precision profile plot with fitted linear regression model with lowest RMSE for lot N001.

4.3 LoQ Results for non-linear models

This analysis fitted multiple non-linear models to the precision profile of within-laboratory variance (SD^2_{WL}) versus mean concentration of the low level samples. SD^2_{WL} used for fitting the models was transformed into CV_{WL} for plotting the precision profile (Figure 2). These models are variance functions commonly applied in clinical literature.

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1
		Page 7 of 24

Table 5 summarizes the evaluated models with their equations, model parameters, Akaike Information Criterion (AIC) and estimated LoQ as the mean concentration at the CV_{WL} goal, sorted from lowest to greatest AIC. Estimated LoQ at CV_{WL} goal (last column, Table 5) as N/A indicates a negative estimated LoQ, that the LoQ cannot be estimated or that the model cannot be estimated (e.g., due to convergence errors or an insufficient number of low level samples). Figure 2 shows the fitted line of the lowest AIC model along with its estimated LoQ.

The multiple fitted models had different numbers of parameters. The proposed model was selected using AIC instead of RMSE since AIC accounts not only for the goodness of fit, but also for model complexity, favoring simpler models with fewer parameters. The lower the AIC, the better, where negative numbers are allowed.

The non-linear precision profile model with the lowest AIC interpolated a LoQ of 0.59 activity.

Table 5. List of the non-linear models fitted to the precision profile for lot N001.

Equation	β_1	β_2	β_3	J	AIC	LoQ at CV_{WL} goal (activity)
$SD_{WL}^2 = \beta_1 \cdot \text{Mean}^l$	0.02111	N/A	N/A	0.8057	- 1709.006	0.59
$SD_{WL}^2 = \beta_1 + \beta_2 \cdot \text{Mean}^l$	9.867e-39	0.02111	N/A	0.8057	- 1709.006	0.59
$SD_{WL}^2 = (\beta_1 + \beta_2 \cdot \text{Mean})^l$	4.297e-17	0.008329	N/A	0.8057	- 1709.006	0.59
$SD_{WL}^2 = (\beta_1 + \beta_2 \cdot \text{Mean})^2$	0.07847	0.05174	N/A	N/A	- 1540.495	0.53
$SD_{WL}^2 = \beta_1 + \beta_2 \cdot \text{Mean}^2$	0.01109	0.004823	N/A	N/A	- 1435.399	0.56
$SD_{WL}^2 = \beta_1$	0.03546	N/A	N/A	N/A	- 1167.802	0.94
$SD_{WL}^2 = \beta_1 + \beta_2 \cdot \text{Mean} + \beta_3 \cdot \text{Mean}^l$	0.01826	- 5.323e+14	5.323e+14	1	- 974.228	22.36*

Abbreviations: N/A, Not Applicable; AIC, Akaike Information Criterion; CV, Coefficient of Variation; WL, Within-Laboratory precision.

*Model extrapolated the LoQ value.

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 8 of 24

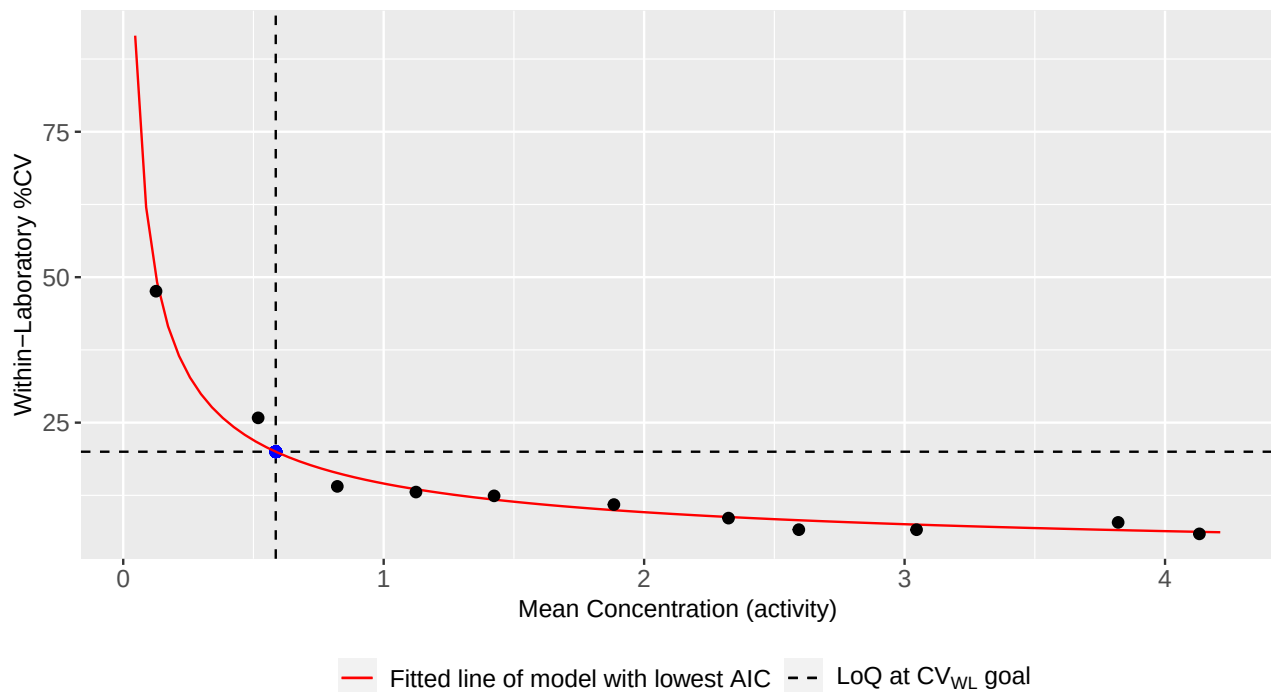


Figure 2. Precision profile plot with fitted non-linear model with lowest AIC for lot N001.

4.4 Results Summary

Claiming a LoQ using a precision profile model depends on selecting an appropriate fit, either by linear or non-linear regression. Visual inspection (Figures 1 and 2) and summary statistics (Tables 4 and 5) help justify model choice. A good model has data points evenly distributed around the fitting line, while gaps or patterns in the data can distort results. When RMSE or AIC values are similar, choose the simpler model with fewer parameters. For linear regression, an R² much lower than 0.9–1 indicates poor fit.

4.4.1 Simple linear regression with variables transformations

The linear precision profile model with lowest RMSE interpolated a LoQ of 0.58 activity. If the fitting looks good, the estimated LoQ might be claimed as the assay’s LoQ. If the fitting does not look good, it is not recommended to claim the interpolated LoQ as the assay’s LoQ. Revisit Table 4 in case other fitted models interpolated the LoQ while keeping a good fitting.

Alternatively, sort low-level sample results from lowest to highest and select the lowest result that meets the precision goal, ensuring all higher results also satisfy this goal. As a worst-case scenario, select the highest low-level sample result that meets the precision goal.

	Statistical Data Analysis Report Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1
		Page 9 of 24

4.4.2 Non-linear models

The non-linear precision profile model with lowest AIC interpolated a LoQ of 0.59 activity. If the fitting looks good, the estimated LoQ might be claimed as the assay's LoQ. If the fitting does not look good, it is not recommended to claim the interpolated LoQ as the assay's LoQ. Revisit Table 5 in case other fitted models interpolated the LoQ while keeping a good fitting.

Alternatively, sort low-level sample results from lowest to highest and select the lowest result that meets the precision goal, ensuring all higher results also satisfy this goal. As a worst-case scenario, select the highest low-level sample result that meets the precision goal.

5 Software

This report has been generated using the following software:

- ShinyMirage version Development
- R version 4.2.2 (2022-10-31 ucrt)
- *openxlsx* package version 4.2.8.1
- *readxl* package version 1.5.0
- *lubridate* package version 1.9.2
- *forcats* package version 1.0.0
- *readr* package version 2.1.4
- *tibble* package version 3.2.1
- *tidyverse* package version 2.0.0
- *kableExtra* package version 1.3.4
- *knitr* package version 1.44
- *purrr* package version 1.0.2
- *VFP* package version 1.4.4
- *gnm* package version 1.1.5
- *tidyr* package version 1.3.0
- *gridExtra* package version 2.3
- *VCA* package version 1.5.1
- *lmtest* package version 0.9.40

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1
		Page 10 of 24

- zoo package version 1.8.11
- reshape2 package version 1.4.4
- reshape package version 0.8.9
- dplyr package version 1.1.3
- zip package version 2.3.0
- qqplotr package version 0.0.6
- ggplot2 package version 3.4.3
- stringr package version 1.5.0
- shinythemes package version 1.2.0
- shiny package version 1.7.5
- stats package version 4.2.2
- graphics package version 4.2.2
- grDevices package version 4.2.2
- datasets package version 4.2.2
- utils package version 4.2.2
- methods package version 4.2.2
- base package version 4.2.2

6 Operator and Data Input Reviewer

Operators:

- Joey (6/25/2025)

Data Input Reviewers:

- Joe (6/25/2025)
- Joseph (6/25/2025)

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 11 of 24

7 Summary Raw Data

Data entry was performed by two different data input reviewers via two different Excel files. Previous to the data analysis, application automatically verified that data in both Excel files were identical.

Table 6. Summary of the Raw Data for lot N001

Sample	Day	Run	Replicate	Result (activity)
DP	1	1	Rep1	0.10
DP	1	1	Rep2	0.10
DP	1	1	Rep3	0.10
DP	1	1	Rep4	0.20
DP	1	1	Rep5	0.10
DP	1	1	Rep6	0.10
DP	1	1	Rep7	0.10
DP	1	1	Rep8	0.10
DP	1	1	Rep9	0.10
DP	1	1	Rep10	0.20
DP	2	1	Rep1	0.10
DP	2	1	Rep2	0.20
DP	2	1	Rep3	0.10
DP	2	1	Rep4	0.00
DP	2	1	Rep5	0.10
DP	2	1	Rep6	0.10
DP	2	1	Rep7	0.10
DP	2	1	Rep8	0.10
DP	2	1	Rep9	0.10
DP	2	1	Rep10	0.10
DP	3	1	Rep1	0.10
DP	3	1	Rep2	0.20
DP	3	1	Rep3	0.20
DP	3	1	Rep4	0.20
DP	3	1	Rep5	0.20
DP	3	1	Rep6	0.00
DP	3	1	Rep7	0.20
DP	3	1	Rep8	0.10
DP	3	1	Rep9	0.00
DP	3	1	Rep10	0.20
DP	4	1	Rep1	0.20
DP	4	1	Rep2	0.20
DP	4	1	Rep3	0.20
DP	4	1	Rep4	0.20
DP	4	1	Rep5	0.10
DP	4	1	Rep6	0.10

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 12 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
DP	4	1	Rep7	0.10
DP	4	1	Rep8	0.10
DP	4	1	Rep9	0.10
DP	4	1	Rep10	0.10
DP	5	1	Rep1	0.20
DP	5	1	Rep2	0.20
DP	5	1	Rep3	0.20
DP	5	1	Rep4	0.20
DP	5	1	Rep5	0.10
DP	5	1	Rep6	0.00
DP	5	1	Rep7	0.10
DP	5	1	Rep8	0.10
DP	5	1	Rep9	0.10
DP	5	1	Rep10	0.10
A	1	1	Rep1	0.50
A	1	1	Rep2	0.50
A	1	1	Rep3	0.50
A	1	1	Rep4	0.50
A	1	1	Rep5	0.50
A	1	1	Rep6	0.80
A	1	1	Rep7	0.50
A	1	1	Rep8	0.50
A	1	1	Rep9	0.50
A	1	1	Rep10	0.50
A	2	1	Rep1	0.60
A	2	1	Rep2	0.50
A	2	1	Rep3	0.40
A	2	1	Rep4	0.60
A	2	1	Rep5	0.50
A	2	1	Rep6	0.40
A	2	1	Rep7	0.50
A	2	1	Rep8	0.40
A	2	1	Rep9	0.40
A	2	1	Rep10	0.50
A	3	1	Rep1	0.50
A	3	1	Rep2	0.50
A	3	1	Rep3	0.50
A	3	1	Rep4	0.70
A	3	1	Rep5	0.80
A	3	1	Rep6	0.30
A	3	1	Rep7	0.30
A	3	1	Rep8	0.50

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 13 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
A	3	1	Rep9	0.50
A	3	1	Rep10	0.70
A	4	1	Rep1	0.70
A	4	1	Rep2	0.60
A	4	1	Rep3	0.70
A	4	1	Rep4	0.70
A	4	1	Rep5	0.80
A	4	1	Rep6	0.50
A	4	1	Rep7	0.40
A	4	1	Rep8	0.50
A	4	1	Rep9	0.60
A	4	1	Rep10	0.60
A	5	1	Rep1	0.60
A	5	1	Rep2	0.50
A	5	1	Rep3	0.50
A	5	1	Rep4	0.70
A	5	1	Rep5	0.50
A	5	1	Rep6	0.40
A	5	1	Rep7	0.30
A	5	1	Rep8	0.30
A	5	1	Rep9	0.30
A	5	1	Rep10	0.30
B	1	1	Rep1	0.70
B	1	1	Rep2	0.90
B	1	1	Rep3	0.70
B	1	1	Rep4	0.90
B	1	1	Rep5	0.80
B	1	1	Rep6	0.80
B	1	1	Rep7	0.90
B	1	1	Rep8	0.70
B	1	1	Rep9	0.90
B	1	1	Rep10	0.90
B	2	1	Rep1	0.90
B	2	1	Rep2	0.80
B	2	1	Rep3	0.70
B	2	1	Rep4	0.80
B	2	1	Rep5	0.80
B	2	1	Rep6	0.50
B	2	1	Rep7	0.80
B	2	1	Rep8	0.80
B	2	1	Rep9	0.70
B	2	1	Rep10	0.80

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 14 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
B	3	1	Rep1	0.80
B	3	1	Rep2	0.80
B	3	1	Rep3	0.90
B	3	1	Rep4	1.00
B	3	1	Rep5	0.80
B	3	1	Rep6	0.80
B	3	1	Rep7	1.10
B	3	1	Rep8	0.60
B	3	1	Rep9	0.90
B	3	1	Rep10	0.90
B	4	1	Rep1	0.90
B	4	1	Rep2	0.90
B	4	1	Rep3	1.00
B	4	1	Rep4	1.00
B	4	1	Rep5	1.00
B	4	1	Rep6	0.80
B	4	1	Rep7	0.90
B	4	1	Rep8	0.80
B	4	1	Rep9	0.70
B	4	1	Rep10	0.70
B	5	1	Rep1	1.00
B	5	1	Rep2	0.90
B	5	1	Rep3	0.90
B	5	1	Rep4	0.80
B	5	1	Rep5	0.80
B	5	1	Rep6	0.70
B	5	1	Rep7	0.70
B	5	1	Rep8	0.70
B	5	1	Rep9	0.70
B	5	1	Rep10	0.80
C	1	1	Rep1	1.20
C	1	1	Rep2	1.20
C	1	1	Rep3	0.90
C	1	1	Rep4	1.40
C	1	1	Rep5	1.30
C	1	1	Rep6	0.90
C	1	1	Rep7	1.10
C	1	1	Rep8	1.20
C	1	1	Rep9	1.20
C	1	1	Rep10	1.10
C	2	1	Rep1	1.00
C	2	1	Rep2	1.20

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 15 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
C	2	1	Rep3	1.00
C	2	1	Rep4	1.10
C	2	1	Rep5	1.10
C	2	1	Rep6	1.10
C	2	1	Rep7	1.00
C	2	1	Rep8	1.00
C	2	1	Rep9	0.90
C	2	1	Rep10	1.10
C	3	1	Rep1	1.00
C	3	1	Rep2	1.20
C	3	1	Rep3	1.20
C	3	1	Rep4	1.20
C	3	1	Rep5	1.00
C	3	1	Rep6	1.10
C	3	1	Rep7	1.30
C	3	1	Rep8	1.30
C	3	1	Rep9	1.00
C	3	1	Rep10	1.10
C	4	1	Rep1	1.30
C	4	1	Rep2	1.20
C	4	1	Rep3	1.30
C	4	1	Rep4	1.40
C	4	1	Rep5	1.30
C	4	1	Rep6	1.10
C	4	1	Rep7	1.10
C	4	1	Rep8	1.30
C	4	1	Rep9	1.20
C	4	1	Rep10	1.20
C	5	1	Rep1	1.10
C	5	1	Rep2	1.10
C	5	1	Rep3	1.10
C	5	1	Rep4	1.20
C	5	1	Rep5	1.30
C	5	1	Rep6	0.90
C	5	1	Rep7	1.00
C	5	1	Rep8	1.00
C	5	1	Rep9	0.70
C	5	1	Rep10	1.00
D	1	1	Rep1	1.50
D	1	1	Rep2	1.50
D	1	1	Rep3	1.50
D	1	1	Rep4	1.40

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 16 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
D	1	1	Rep5	1.40
D	1	1	Rep6	1.20
D	1	1	Rep7	1.60
D	1	1	Rep8	1.40
D	1	1	Rep9	1.30
D	1	1	Rep10	1.50
D	2	1	Rep1	1.30
D	2	1	Rep2	1.60
D	2	1	Rep3	1.40
D	2	1	Rep4	1.40
D	2	1	Rep5	1.30
D	2	1	Rep6	1.60
D	2	1	Rep7	1.40
D	2	1	Rep8	1.10
D	2	1	Rep9	1.10
D	2	1	Rep10	1.70
D	3	1	Rep1	1.50
D	3	1	Rep2	1.20
D	3	1	Rep3	1.20
D	3	1	Rep4	1.90
D	3	1	Rep5	1.20
D	3	1	Rep6	1.40
D	3	1	Rep7	1.50
D	3	1	Rep8	1.20
D	3	1	Rep9	1.80
D	3	1	Rep10	1.30
D	4	1	Rep1	1.60
D	4	1	Rep2	1.50
D	4	1	Rep3	1.40
D	4	1	Rep4	1.40
D	4	1	Rep5	1.50
D	4	1	Rep6	1.60
D	4	1	Rep7	1.40
D	4	1	Rep8	1.40
D	4	1	Rep9	1.70
D	4	1	Rep10	1.60
D	5	1	Rep1	1.40
D	5	1	Rep2	1.60
D	5	1	Rep3	1.60
D	5	1	Rep4	1.20
D	5	1	Rep5	1.50
D	5	1	Rep6	1.20

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 17 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
D	5	1	Rep7	1.30
D	5	1	Rep8	1.30
D	5	1	Rep9	1.40
D	5	1	Rep10	1.20
E	1	1	Rep1	1.80
E	1	1	Rep2	1.90
E	1	1	Rep3	1.80
E	1	1	Rep4	2.00
E	1	1	Rep5	2.00
E	1	1	Rep6	2.10
E	1	1	Rep7	1.50
E	1	1	Rep8	1.80
E	1	1	Rep9	1.90
E	1	1	Rep10	1.60
E	2	1	Rep1	1.70
E	2	1	Rep2	1.90
E	2	1	Rep3	1.70
E	2	1	Rep4	1.90
E	2	1	Rep5	1.80
E	2	1	Rep6	1.30
E	2	1	Rep7	1.80
E	2	1	Rep8	1.90
E	2	1	Rep9	1.90
E	2	1	Rep10	2.00
E	3	1	Rep1	1.90
E	3	1	Rep2	1.80
E	3	1	Rep3	2.10
E	3	1	Rep4	1.70
E	3	1	Rep5	2.20
E	3	1	Rep6	2.10
E	3	1	Rep7	1.60
E	3	1	Rep8	1.80
E	3	1	Rep9	2.00
E	3	1	Rep10	2.10
E	4	1	Rep1	1.80
E	4	1	Rep2	1.80
E	4	1	Rep3	2.00
E	4	1	Rep4	2.00
E	4	1	Rep5	1.80
E	4	1	Rep6	2.10
E	4	1	Rep7	2.10
E	4	1	Rep8	2.20

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 18 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
E	4	1	Rep9	2.00
E	4	1	Rep10	2.20
E	5	1	Rep1	2.10
E	5	1	Rep2	2.10
E	5	1	Rep3	2.00
E	5	1	Rep4	2.10
E	5	1	Rep5	1.50
E	5	1	Rep6	2.00
E	5	1	Rep7	1.60
E	5	1	Rep8	1.90
E	5	1	Rep9	1.80
E	5	1	Rep10	1.50
F	1	1	Rep1	2.20
F	1	1	Rep2	2.80
F	1	1	Rep3	2.20
F	1	1	Rep4	2.40
F	1	1	Rep5	2.20
F	1	1	Rep6	2.00
F	1	1	Rep7	2.50
F	1	1	Rep8	2.00
F	1	1	Rep9	2.50
F	1	1	Rep10	2.80
F	2	1	Rep1	2.40
F	2	1	Rep2	2.50
F	2	1	Rep3	2.30
F	2	1	Rep4	2.20
F	2	1	Rep5	2.10
F	2	1	Rep6	2.20
F	2	1	Rep7	2.50
F	2	1	Rep8	2.00
F	2	1	Rep9	2.30
F	2	1	Rep10	2.20
F	3	1	Rep1	2.30
F	3	1	Rep2	2.50
F	3	1	Rep3	2.30
F	3	1	Rep4	2.60
F	3	1	Rep5	2.30
F	3	1	Rep6	2.20
F	3	1	Rep7	2.40
F	3	1	Rep8	2.20
F	3	1	Rep9	2.50
F	3	1	Rep10	2.60

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 19 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
F	4	1	Rep1	2.30
F	4	1	Rep2	2.60
F	4	1	Rep3	2.10
F	4	1	Rep4	2.30
F	4	1	Rep5	2.30
F	4	1	Rep6	2.40
F	4	1	Rep7	2.40
F	4	1	Rep8	2.50
F	4	1	Rep9	2.50
F	4	1	Rep10	1.90
F	5	1	Rep1	2.30
F	5	1	Rep2	2.30
F	5	1	Rep3	2.10
F	5	1	Rep4	2.30
F	5	1	Rep5	2.30
F	5	1	Rep6	2.30
F	5	1	Rep7	2.40
F	5	1	Rep8	1.90
F	5	1	Rep9	2.30
F	5	1	Rep10	2.50
G	1	1	Rep1	2.50
G	1	1	Rep2	2.70
G	1	1	Rep3	2.50
G	1	1	Rep4	2.80
G	1	1	Rep5	2.70
G	1	1	Rep6	2.60
G	1	1	Rep7	2.50
G	1	1	Rep8	2.60
G	1	1	Rep9	2.20
G	1	1	Rep10	2.60
G	2	1	Rep1	2.60
G	2	1	Rep2	2.50
G	2	1	Rep3	2.40
G	2	1	Rep4	2.50
G	2	1	Rep5	2.50
G	2	1	Rep6	2.50
G	2	1	Rep7	2.50
G	2	1	Rep8	2.40
G	2	1	Rep9	2.60
G	2	1	Rep10	2.40
G	3	1	Rep1	2.60
G	3	1	Rep2	2.60

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 20 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
G	3	1	Rep3	2.40
G	3	1	Rep4	2.90
G	3	1	Rep5	2.90
G	3	1	Rep6	2.40
G	3	1	Rep7	2.40
G	3	1	Rep8	2.90
G	3	1	Rep9	2.90
G	3	1	Rep10	2.80
G	4	1	Rep1	2.60
G	4	1	Rep2	2.60
G	4	1	Rep3	2.60
G	4	1	Rep4	2.80
G	4	1	Rep5	2.70
G	4	1	Rep6	2.70
G	4	1	Rep7	3.00
G	4	1	Rep8	2.60
G	4	1	Rep9	2.50
G	4	1	Rep10	2.50
G	5	1	Rep1	2.50
G	5	1	Rep2	2.70
G	5	1	Rep3	2.50
G	5	1	Rep4	2.70
G	5	1	Rep5	2.90
G	5	1	Rep6	2.40
G	5	1	Rep7	2.50
G	5	1	Rep8	2.50
G	5	1	Rep9	2.40
G	5	1	Rep10	2.60
H	1	1	Rep1	2.90
H	1	1	Rep2	3.10
H	1	1	Rep3	3.00
H	1	1	Rep4	3.30
H	1	1	Rep5	3.40
H	1	1	Rep6	3.30
H	1	1	Rep7	2.80
H	1	1	Rep8	3.00
H	1	1	Rep9	3.00
H	1	1	Rep10	3.20
H	2	1	Rep1	2.70
H	2	1	Rep2	2.90
H	2	1	Rep3	2.70
H	2	1	Rep4	2.90

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 21 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
H	2	1	Rep5	2.90
H	2	1	Rep6	3.10
H	2	1	Rep7	3.40
H	2	1	Rep8	2.90
H	2	1	Rep9	2.90
H	2	1	Rep10	3.10
H	3	1	Rep1	3.20
H	3	1	Rep2	3.10
H	3	1	Rep3	3.10
H	3	1	Rep4	2.90
H	3	1	Rep5	3.30
H	3	1	Rep6	2.80
H	3	1	Rep7	3.10
H	3	1	Rep8	3.00
H	3	1	Rep9	3.00
H	3	1	Rep10	3.30
H	4	1	Rep1	3.00
H	4	1	Rep2	3.30
H	4	1	Rep3	3.30
H	4	1	Rep4	3.00
H	4	1	Rep5	3.10
H	4	1	Rep6	3.00
H	4	1	Rep7	3.20
H	4	1	Rep8	3.00
H	4	1	Rep9	3.30
H	4	1	Rep10	3.20
H	5	1	Rep1	3.40
H	5	1	Rep2	2.70
H	5	1	Rep3	3.10
H	5	1	Rep4	3.00
H	5	1	Rep5	3.10
H	5	1	Rep6	3.00
H	5	1	Rep7	2.90
H	5	1	Rep8	3.10
H	5	1	Rep9	2.60
H	5	1	Rep10	2.70
I	1	1	Rep1	3.90
I	1	1	Rep2	3.50
I	1	1	Rep3	3.70
I	1	1	Rep4	3.80
I	1	1	Rep5	4.50
I	1	1	Rep6	3.00

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 22 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
I	1	1	Rep7	3.80
I	1	1	Rep8	3.40
I	1	1	Rep9	3.70
I	1	1	Rep10	4.30
I	2	1	Rep1	3.70
I	2	1	Rep2	3.70
I	2	1	Rep3	3.60
I	2	1	Rep4	3.70
I	2	1	Rep5	3.70
I	2	1	Rep6	3.90
I	2	1	Rep7	3.50
I	2	1	Rep8	4.00
I	2	1	Rep9	3.40
I	2	1	Rep10	4.10
I	3	1	Rep1	3.70
I	3	1	Rep2	4.10
I	3	1	Rep3	3.70
I	3	1	Rep4	3.30
I	3	1	Rep5	3.90
I	3	1	Rep6	3.60
I	3	1	Rep7	4.00
I	3	1	Rep8	4.20
I	3	1	Rep9	3.70
I	3	1	Rep10	4.00
I	4	1	Rep1	4.20
I	4	1	Rep2	3.80
I	4	1	Rep3	3.70
I	4	1	Rep4	4.00
I	4	1	Rep5	4.00
I	4	1	Rep6	4.20
I	4	1	Rep7	4.10
I	4	1	Rep8	4.10
I	4	1	Rep9	4.40
I	4	1	Rep10	4.00
I	5	1	Rep1	3.40
I	5	1	Rep2	3.70
I	5	1	Rep3	3.40
I	5	1	Rep4	3.90
I	5	1	Rep5	3.80
I	5	1	Rep6	4.10
I	5	1	Rep7	3.90
I	5	1	Rep8	3.60

werfen	Statistical Data Analysis Report	
	Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1 Page 23 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
I	5	1	Rep9	3.60
I	5	1	Rep10	4.00
J	1	1	Rep1	4.60
J	1	1	Rep2	4.20
J	1	1	Rep3	4.30
J	1	1	Rep4	4.40
J	1	1	Rep5	4.10
J	1	1	Rep6	4.20
J	1	1	Rep7	3.90
J	1	1	Rep8	4.20
J	1	1	Rep9	4.50
J	1	1	Rep10	4.40
J	2	1	Rep1	3.70
J	2	1	Rep2	4.10
J	2	1	Rep3	3.90
J	2	1	Rep4	4.10
J	2	1	Rep5	3.90
J	2	1	Rep6	4.30
J	2	1	Rep7	4.00
J	2	1	Rep8	4.10
J	2	1	Rep9	4.50
J	2	1	Rep10	4.00
J	3	1	Rep1	3.90
J	3	1	Rep2	4.10
J	3	1	Rep3	4.00
J	3	1	Rep4	3.80
J	3	1	Rep5	4.20
J	3	1	Rep6	4.10
J	3	1	Rep7	4.30
J	3	1	Rep8	4.00
J	3	1	Rep9	3.90
J	3	1	Rep10	4.50
J	4	1	Rep1	4.10
J	4	1	Rep2	3.90
J	4	1	Rep3	4.10
J	4	1	Rep4	4.40
J	4	1	Rep5	4.10
J	4	1	Rep6	4.30
J	4	1	Rep7	4.70
J	4	1	Rep8	3.60
J	4	1	Rep9	4.60
J	4	1	Rep10	3.80

werfen	Statistical Data Analysis Report Study: LoQ	
	Assay name: FVIII_Low_Test Platform: ModTOP Project code: BEAHCA34 Time Stamp: 2026-05-26, 10:40:33 (UTC)	Edition 1
		Page 24 of 24

Table 6. Summary of the Raw Data for lot N001 (*continued*)

Sample	Day	Run	Replicate	Result (activity)
J	5	1	Rep1	4.00
J	5	1	Rep2	4.10
J	5	1	Rep3	4.10
J	5	1	Rep4	4.10
J	5	1	Rep5	4.00
J	5	1	Rep6	4.30
J	5	1	Rep7	4.40
J	5	1	Rep8	3.90
J	5	1	Rep9	3.90
J	5	1	Rep10	4.00